

Morrison

STEERABLE NEEDLE®

Your first choice for
deep-seated targets



Indication for use:

Percutaneous injection and aspiration procedures and tissue sampling in kidney, liver, muscle and subcutaneous tissue.



AprioMed

- **Steer distal end of needle with live feedback under image guidance**
- **Facilitates precise needle placement**
- **Curve around sensitive structures or other obstacles**



Deep-seated target access

Morrison Steerable Needle™ used for access to a deep-seated target in the hip joint.

USE MORRISON STEERABLE NEEDLE® TO:


- Make major adjustments around bones, organs or neurovascular structures
- Make minor adjustments near target
- Make multiple adjustments as needle is advanced

SEAMLESS ADJUSTMENT

Gives precise control

FUNCTION

When lever is turned, needle curves towards lever

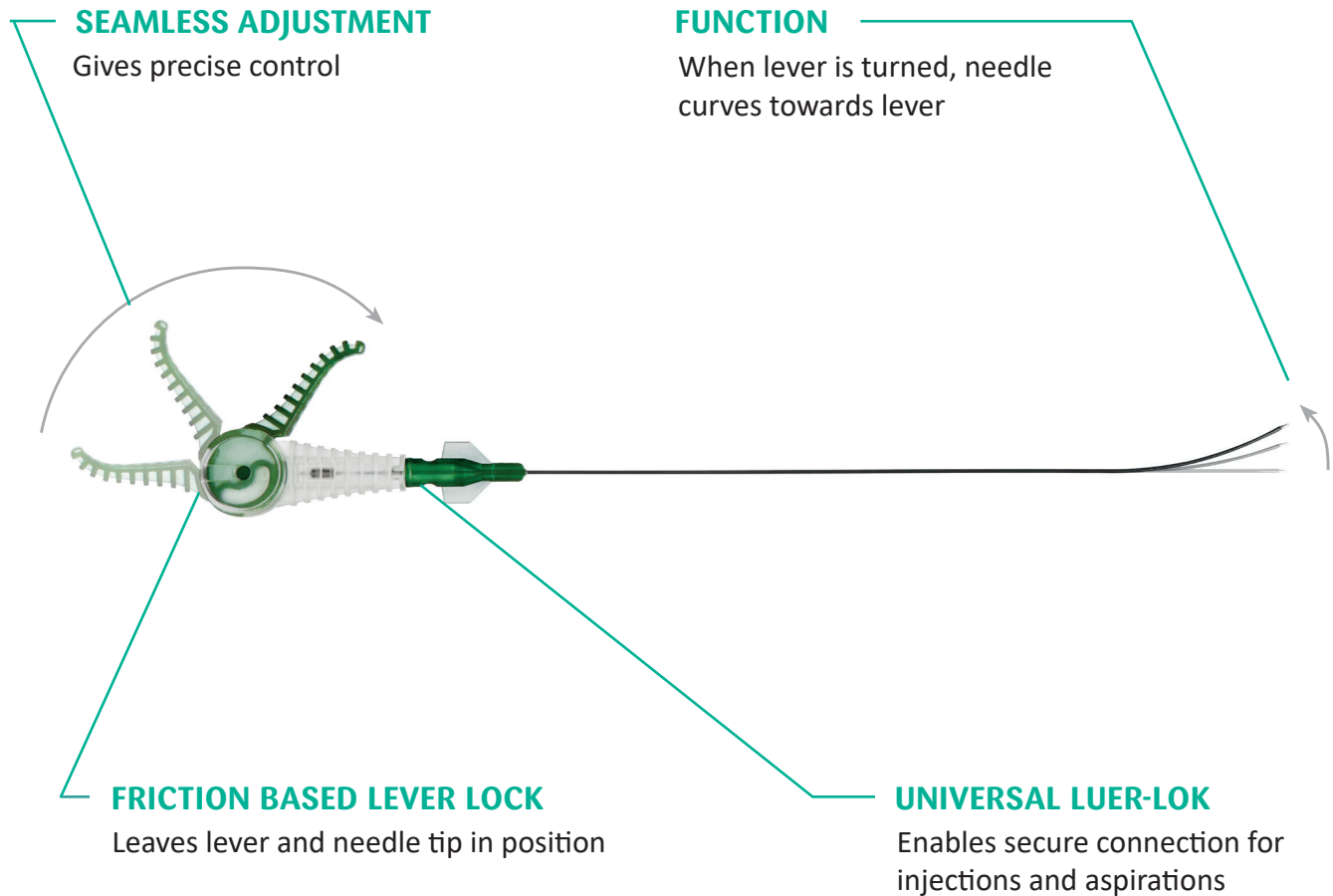

FRICITION BASED LEVER LOCK

Leaves lever and needle tip in position

UNIVERSAL LUER-LOK


Enables secure connection for injections and aspirations

PENCIL POINT STYLET TIP

Facilitates straight and smooth needle advancement

30° BACK BEVEL NEEDLE TIP

Provides perfect cut through tissue

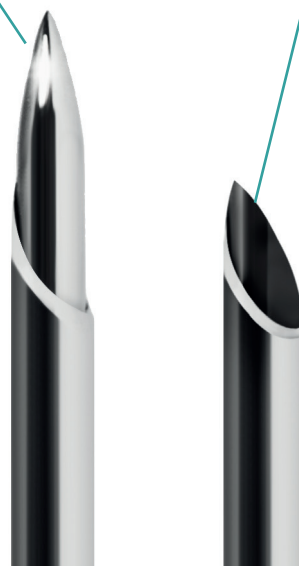


TECHNICAL SPECIFICATIONS

Morrison Steerable Needle®, 21G


Package contents:

- Steerable Stylet
- Needle
 - Inner diameter 0.68 mm
 - Outer diameter 0.82 mm (21G)
 - Length 167 mm

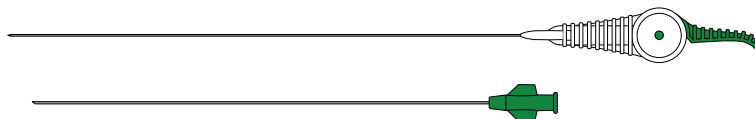


PERFORMANCE SPECIFICATIONS

- Length, upright position
- Deflectable tip length
- Maximum tip deflection
- Curve radius
- Seamless adjustment

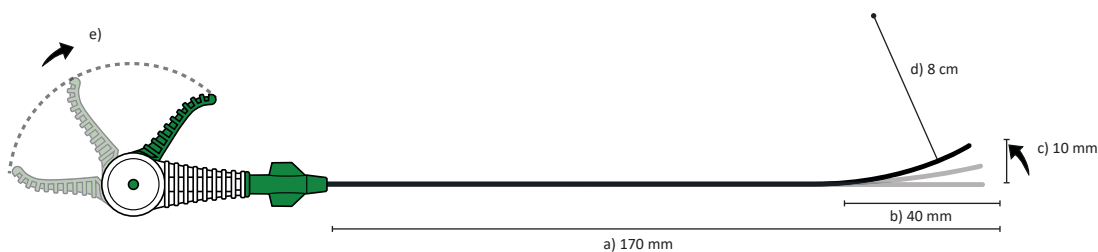

ORDERING INFORMATION

Product number	Product description	Minimum Order Quantity
13550	Morrison Steerable Needle® 21G 17 cm	10 units

 Indication for use: Percutaneous injection and aspiration procedures and tissue sampling in kidney, liver, muscle and subcutaneous tissue. Before using Morrison Steerable Needle® read the instructions for use which accompany the product for indications, contraindications, warnings and precautions. Morrison Steerable Needle® is a trademark of AprioMed AB. Patents pending.

Available for sale in the USA only.

HEAD OFFICE

AprioMed AB
Viridings Allé 28
SE-754 50 Uppsala, Sweden
Phone: +46 18 430 14 40
Free Fax (Int'l): +46 770 339972

US SALES OFFICE

AprioMed, Inc.
45 South Main St Unit 2
Derry, NH 03038, USA
Phone: +1 603 421 0875
Fax (US toll free): 866 239 179



info@apriomed.com / www.apriomed.com

APRIOMED develops, manufactures, markets and sells innovative medical devices and related services within the field of interventional radiology. We aim to deliver, in close collaboration with healthcare practitioners, innovative quality tools to achieve the most favourable solutions for radiologists worldwide. AprioMed is headquartered in Uppsala, Sweden, with a sales subsidiary in Londonderry, NH, USA. Our products are sold in more than 40 countries by AprioMed (direct sales in US, UK, Ireland and Sweden) and our local distributors.