

Breast Lesion Localization / Kopans

Kopans

Breast Lesion Localization Needle

Used for preoperative marking of nonpalpable breast lesions

- The spring-hook wire acts as a marker to **indicate the location** of the lesion
- The reference burnish mark, when visible just outside the needle hub on the spring-hook wire, **provides visual assurance** that the hook wire is within the needle tip during needle manipulation
- The reinforced portion of the modified spring-hook wire serves as a **visual and palpable landmark** that helps the surgeon gauge the distance to the lesion during dissection and **reduces the risk** of severing the spring-hook wire
- The EchoTip® models allow for **improved visibility** under ultrasound

Supplied sterile in peel-open packages. Intended for one-time use only.




ORDER INFORMATION

Kopans - Standard Needles

SET ORDER NUMBER	OUTER CANNULA GAUGE • LENGTH	GLOBAL PROD. NO.
DKBL-20-5.0-A	20G • 5 cm	G02586
DKBL-20-7.0-A	20G • 7 cm	G02601
DKBL-20-9.0-A	20G • 9 cm	G02587
DKBL-20-15.0-A	20G • 15 cm	G02850
DKBL-21-5.0-A	21G • 5 cm	G02900
DKBL-21-7.0-A	21G • 7 cm	G03507
DKBL-21-9.0-A	21G • 9 cm	G02554
DKBL-21-15.0-A	21G • 15 cm	G02938

Kopans - EchoTip® Needle

SET ORDER NUMBER	OUTER CANNULA GAUGE • LENGTH	GLOBAL PROD. NO.
DKBL-20-5.0-A-U	20G • 5 cm	G04343
DKBL-20-7.0-A-U	20G • 7 cm	G09984
DKBL-20-9.0-A-U	20G • 9 cm	G04296
DKBL-21-5.0-A-U	21G • 5 cm	G05862
DKBL-21-7.0-A-U	21G • 7 cm	G09042
DKBL-21-9.0-A-U	21G • 9 cm	G07370

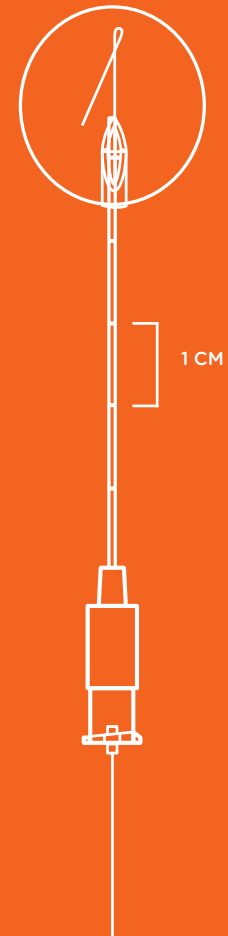

Not all parts numbers shown may be approved for sale in all regulatory jurisdictions. Parts and configurations are subject to change without notification. Please consult with your Account Representative for details.