



Osteo-Site® Products / Cement Injectors

Duro-Ject®

Bone Cement Injector Set

Used for injection of low, medium and high viscosity bone cement into surgical sites.

Set includes:

- 2 rotating Luer Lock connectors proving free rotational movement of the injector

Supplied sterile in peel-open packages. Intended for one-time use only.



FEATURES & BENEFITS

- Injector incorporates an ergonomic design, allowing uncomplicated operation
- Injector is supplied with two 10ml injector barrels that are graduated for easy monitoring of low and medium viscosity cement flow
- Set is supplied with two 29cm precurved connecting tubes and two cement funnels for easy loading
- Precision-ground metal threads improve accuracy and control of cement injection
- Each full rotation of injector handle delivers 0.2 ml of cement


Enhance your **procedural control** with the **ergonomically robust injector** for **low and medium viscosity** cement.

The included 29cm connecting tube **reduces your radiation exposure**. 1 set **extra is included** for a **bipedicular approach**

ORDER INFORMATION

Duro-Ject®

SET ORDER NUMBER	SYRINGE VOLUME	CONNECTING TUBE LENGTH	CONNECTING TUBE VOLUME	GLOBAL PRODUCT NUMBER
DURO-100	10ml	29cm	1.4ml	G32497

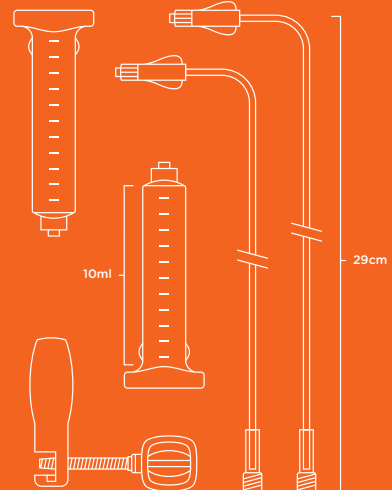

Not all parts numbers shown may be approved for sale in all regulatory jurisdictions. Parts and configurations are subject to change without notification. Please consult with your Account Representative for details.

Osteo-Site® Products / Cement Injectors

Osteo-Force®

Bone Cement Injector Set

Used for injection of low and medium viscosity bone cement into surgical sites.



Supplied sterile in peel-open packages. Intended for one-time use only.

FEATURES & BENEFITS

- Injector has an easy-to-control handle for uncomplicated two-handed operation
- Barrel of syringe holds up to 10 ml of cement
- Syringe incorporates selfpurging technology that evacuates air from the system prior to use
- Supplied with a 15 cm connecting tube
- Each half-rotation of injector handle delivers 0.2 ml of cement


Reduce procedure **time** and **costs** with a lightweight, **easy-to-use** injector.

The selfpurging **injector system** is for **low viscosity** bone cement.

ORDER INFORMATION

Duro-Ject®

SET ORDER NUMBER	SYRINGE VOLUME	CONNECTING TUBE LENGTH	CONNECTING TUBE VOLUME	GLOBAL PRODUCT NUMBER
HPIS-100	10ml	15cm	0.8ml	G13082

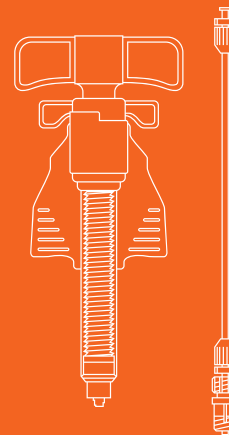

Not all parts numbers shown may be approved for sale in all regulatory jurisdictions. Parts and configurations are subject to change without notification. Please consult with your Account Representative for details.